

UK Office

Everest Biotech Ltd

Cherwell Innovation Centre 77 Heyford Park Upper Heyford Oxfordshire OX25 5HD UK

Enquiries:

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

EB10143 - Goat Anti-FANCE Antibody

Size: 100µg specific antibody in 200µl



Target Protein

Principal Names: FACE, FAE, Fanconi anemia, complementation group E, FANCE

Official Symbol: FANCE

Accession Number(s): NP_068741.1

Human GeneID(s): 2178

Immunogen

Peptide with sequence C-EHKSLESLADGGS, from the internal region of the protein sequence according to NP_068741.1.

Please note the peptide is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested

Peptide ELISA: antibody detection limit dilution 1:32000.


Western blot: Preliminary experiments gave an approx 70kDa band in Human Adipose and Uterus lysates after 0.1μg/ml antibody staining. Please note that currently we cannot find an explanation in the literature for the band we observe given the calculated size of 58.7kDa according to NP_068741.1. The 70kDa band was successfully blocked by incubation with the immunizing peptide. We would appreciate any feedback from people in the field - have any results been reported with other antibodies/lysates? Have any further splice variants/modified forms been reported?

IHC: In paraffin embedded Human Placenta shows nuclear staining in some trophoblasts. Recommended concentration, 3-6µg/ml.

Species Reactivity

Tested: Human

Expected from sequence similarity: Human, Cow

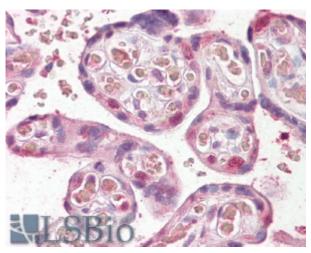

EB10143 (3.8 μ g/ml) staining of paraffin embedded Human Placenta. Steamed antigen retrieval with citrate buffer pH 6, AP-staining.