

#### **UK Office**

#### **Everest Biotech Ltd**

Cherwell Innovation Centre 77 Heyford Park Upper Heyford Oxfordshire OX25 5HD UK

Enquiries:

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

# EB06491 - Goat Anti-FBP17 / FNBP1 Antibody

Size: 100µg specific antibody in 200µl



## **Target Protein**

Principal Names: FNBP1, FBP17, KIAA0554, formin binding protein 1, formin-binding

protein 17, MGC126804

Official Symbol: FNBP1

Accession Number(s): NP\_055848.1

Human GenelD(s): 23048

## **Immunogen**

Peptide with sequence KQLESSKRRFERDC, from the internal region of the protein sequence according to NP\_055848.1.

Please note the peptide is available for sale.

## **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

### **Applications Tested**

Peptide ELISA: antibody detection limit dilution 1:32000.

**Western blot:** Approx 70kDa band observed in Human Testis lysates (calculated MW of 71.3kDa according to NP\_055848.1). Recommended concentration: 0.25-2µg/ml.

**IHC:** Paraffin embedded Human Tonsil and Small Intestine. Recommended concentration: 3.75µg/ml.


### **Species Reactivity**

Tested: Human


Expected from sequence similarity: Human, Mouse, Rat, Dog, Pig, Cow



EB06491 staining (0.25 $\mu$ g/ml) of Human Testis lysate (RIPA buffer, 35 $\mu$ g total protein per lane). Primary incubated for 1 hour. Detected by chemiluminescence.

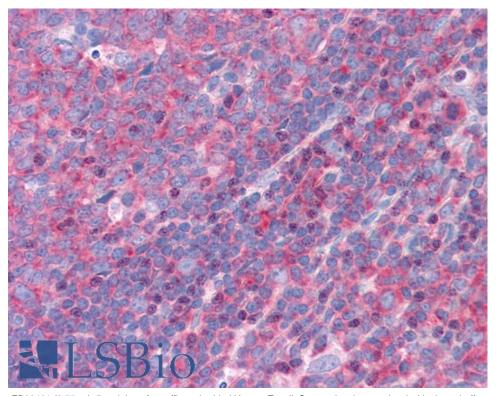

EB06491 (3.75μg/ml) staining of paraffin embedded Human Tonsil. Steamed antigen retrieval with citrate buffer pH 6, AP-staining.

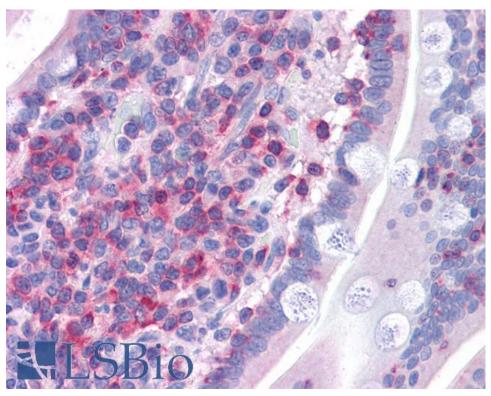

EB06491 (3.75µg/ml) staining of paraffin embedded Human Small Intestine. Steamed antigen retrieval with citrate buffer pH 6, AP-staining.