

#### **UK Office**

#### **Everest Biotech Ltd**

Cherwell Innovation Centre 77 Heyford Park Upper Heyford Oxfordshire OX25 5HD UK

**Enquiries:** 

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

# EB09093 - Goat Anti-GABBR1 Antibody

Size: 100µg specific antibody in 200µl



#### **Target Protein**

**Principal Names:** GABBR1, gamma-aminobutyric acid (GABA) B receptor, 1, FLJ92613, GABAB(1e), GABABR1, GABBR1-3, GPRC3A, dJ271M21.1.1, dJ271M21.1.2, hGB1a, GABA-B receptor, GABAB, subunit 1c, OTTHUMP00000029114,

OTTHUMP0000029115, OTTHUMP0000029116, OTTHUMP00000109099,

gamma-aminobutyric acid (GABA) B receptor 1, seven transmembrane helix receptor

Official Symbol: GABBR1

Accession Number(s): NP\_001461.1; NP\_068703.1; NP\_068704.2; NP\_068705.2

Human GenelD(s): 2550

**Important Comments:** This antibody is expected to recognize reported isoforms a, b and c (NP\_001461.1, NP\_068703.1, NP\_068704.2 resp.).

### **Immunogen**

Peptide with sequence DRLSCDGSRVHLLYK, from the C Terminus of the protein sequence according to NP\_001461.1; NP\_068703.1; NP\_068704.2; NP\_068705.2.

Please note the peptide is available for sale.

## **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

#### **Applications Tested**

Peptide ELISA: antibody detection limit dilution 1:4000.


IHC: In paraffin embedded Human Cortex shows axonal staining. Recommended

concentration: 5-8µg/ml.

#### **Species Reactivity**

Tested: Human

Expected from sequence similarity: Human, Mouse, Rat, Dog, Cow

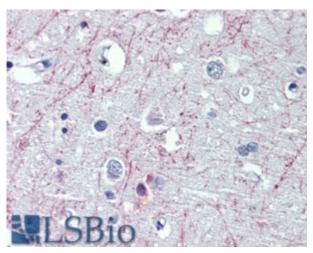

EB09093 (5μg/ml) staining of paraffin embedded Human Cortex. Steamed antigen retrieval with citrate buffer pH 6, AP-staining.