

#### **UK Office**

#### **Everest Biotech Ltd**

Cherwell Innovation Centre 77 Heyford Park Upper Heyford Oxfordshire OX25 5HD UK

Enquiries:

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

# EB06831 - Goat Anti-VPS41 Antibody

Size: 100µg specific antibody in 200µl



## **Target Protein**

**Principal Names:** VPS41, vacuolar protein sorting 41 (yeast), HGNC:12713, HVSP41, vacuolar assembly protein 41, vacuolar protein sorting 41 (yeast homolog), vacuolar protein sorting 41 homolog (S. cerevisiae), HVPS41, hVps41p, vacuolar protein sorting 41

Official Symbol: VPS41

Accession Number(s): NP\_055211.2; NP\_542198.2

Human GeneID(s): 27072

Important Comments: This antibody is expected to recognise both reported isoforms

(NP\_055211.2 and NP\_542198.2)

# Immunogen

Peptide with sequence C-DHLKKDSQNKTLLK, from the internal region of the protein sequence according to NP\_055211.2; NP\_542198.2.

Please note the peptide is available for sale.

## **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin

Aliquot and store at -20°C. Minimize freezing and thawing.

## **Applications Tested**


Peptide ELISA: antibody detection limit dilution 1:16000.

Western blot: Preliminary experiments gave bands at approx 48kDa and 28kDa in human duodenum, kidney and placenta lysates after 1μg/ml antibody staining. Please note that currently we cannot find an explanation in the literature for the bands we observe given the calculated size of 98.5kDa according to NP\_055211.1 and 100kDa according to NP\_542198.1. Both detected bands were successfully blocked by incubation with the immunizing peptide (and BLAST results with the immunizing peptide sequence did not identify any other proteins to explain the additional bands). We would appreciate any feedback from people in the field - have any results been reported with other antibodies/lysates? Have any further splice variants/modified forms been reported IHC: In paraffin embedded Human Kidney shows staining of collective tubules. Recommended concentration, 4-6μg/ml.

### **Species Reactivity**

Tested: Human

Expected from sequence similarity: Human, Mouse, Rat, Dog, Cow

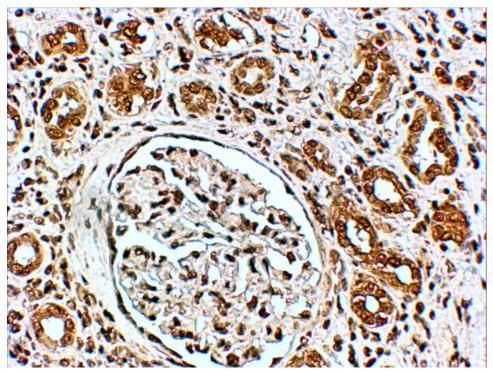

EB06831 (4µg/ml) staining of paraffin embedded Human Kidney. Steamed antigen retrieval with citrate buffer pH 6, HRP-staining. Similar results were obtained after antigen retrieval at pH9.